



medifa H!Light

LED operating light system



we care.

medifa H!Light

- Innovative lighting technology combined with high functionality
- User-friendly control concept
- Operating theatre adapted colour rendering
- Large, homogeneous surface illumination
- Easy adjustability
- One universal handle holder for all cameras (4K and HD) and all handles
- Sterile control at the multifunction handle
- Lightweight construction (aluminium housing)
- Optional: 360° monitor arm with stop-free signal transmission
- Ideal cleaning and disinfection conditions
- Sustainable energy efficiency with minimal warming
- Long lifespan
- Optional with HD or 4K wireless camera



medifa H!Light - illuminating advantages



Light field size

Light field sizes from 18 - 32 cm.

For the distances of 0.8 to 1.2 metres frequently required in everyday operating theatres




LED Technology

The newest medifa LED technology for a homogeneous mixture of cold and warm white LEDs behind the lenses.

The mixed light emerges as white light and is distributed evenly over the wound area.

This prevents the creation of colour shadows by the surgeon's head or hands in the light channel of the operating light.

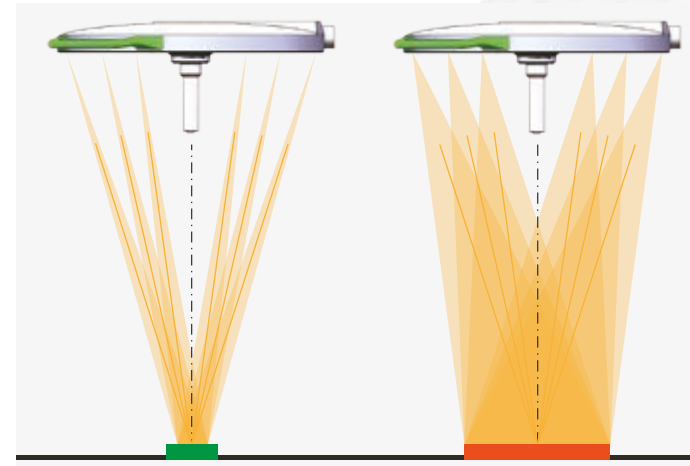

Focus

Electronic focus through the medifa dual lens system, consisting of 72 versatile lenses in 2 sizes.

No mechanical adjustment of the optical components necessary.



Universal handle holder

The universal handle holder offers unique flexibility.

All handle adapters fit into the same handle holder.


Tool-free change possible.

medifa H!Light - illuminating advantages Models



Colour temperature

The right colour temperature for every type of tissue and wound structure.
From 3,000 K to 5,500 K.

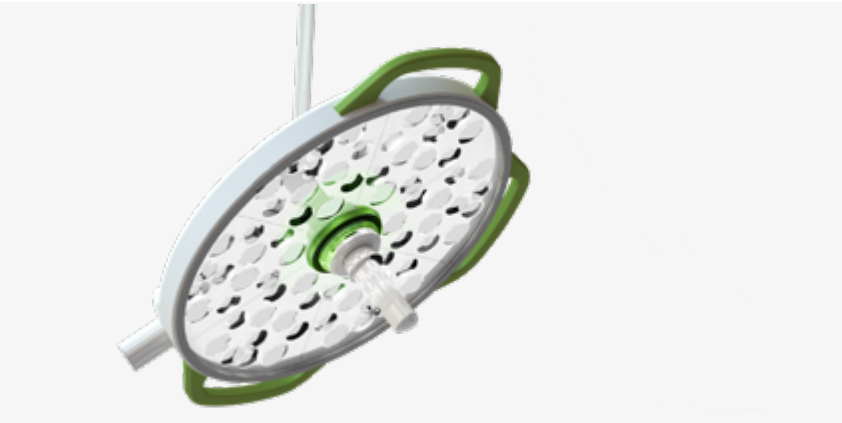

Endoscopy mode

The inner area is illuminated in green for higher contrast on the monitor.
Extremely user-friendly for minimally invasive procedures.



medifa H!Light
Art-No. 100200



medifa H!Light plus
Art-No. 100400

	Colour temperature 3,000 K -5,500 K		
	Control element on light head		
	Sterile control at the multifunction handle		
	Active shadow management		
	Camera		
	Endoscopy mode		
	Memory-function		

included optional




User-friendly control concept

Operating light

Touch control element

on the lamp head and optionally as wall control units surface-mounted (Art. no. 105600), flush-mounted (Art. no. 105650) or panel-mounted (Art. no. 105660)


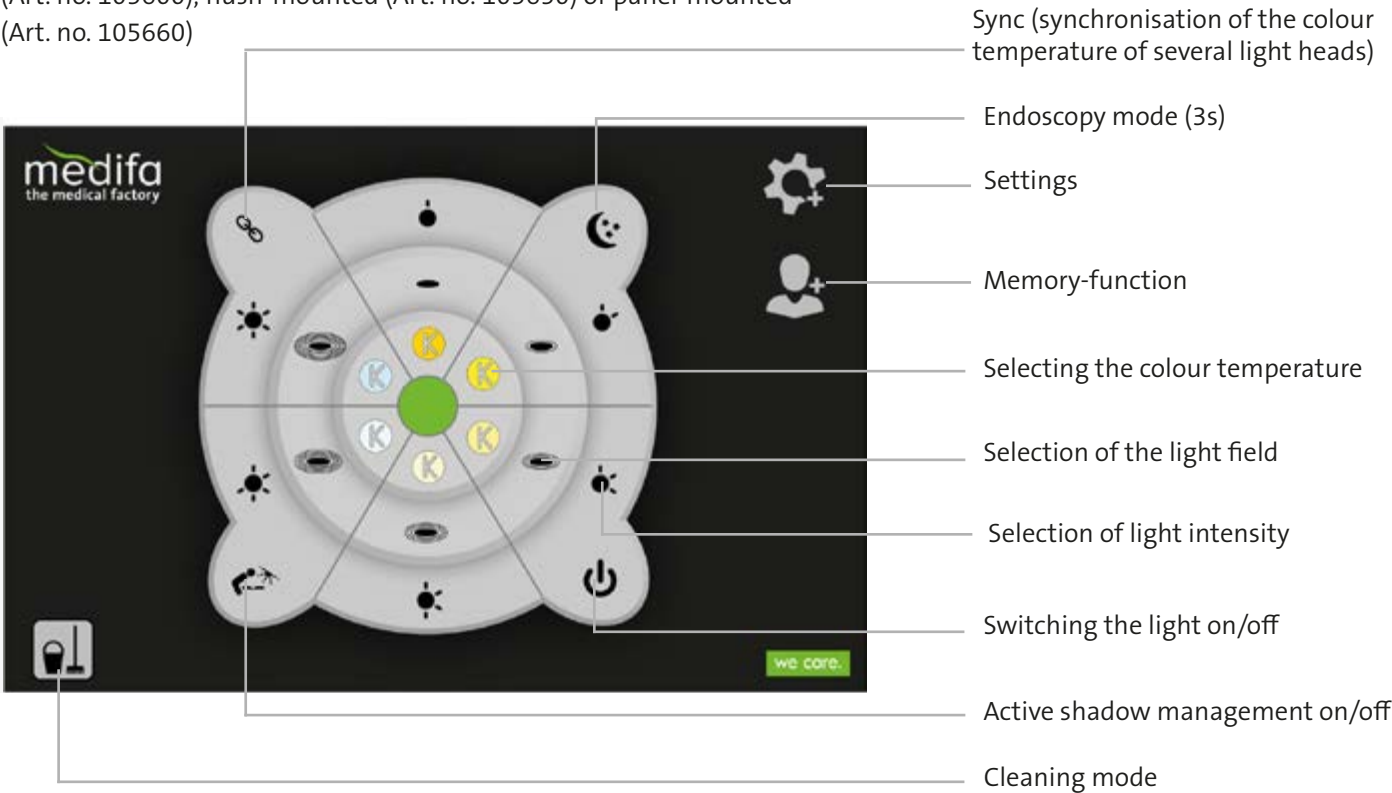


Option: Sterile control at the multifunction handle

Intensity, light field size and colour temperature can be controlled on the sterile handle.
The functions are selected by pressing the selection button at the bottom of the handle.
The selected function is adjusted by turning the sterilisable handle.

Wireless-camera

Remote control

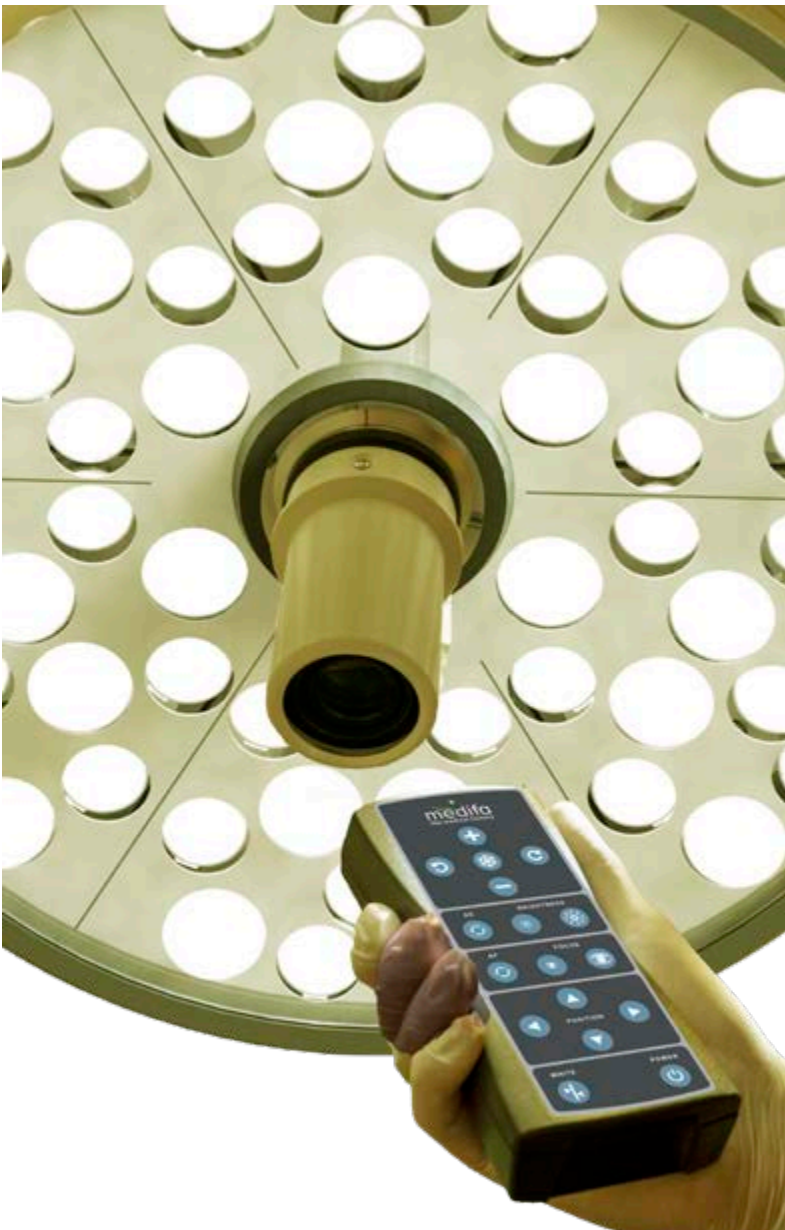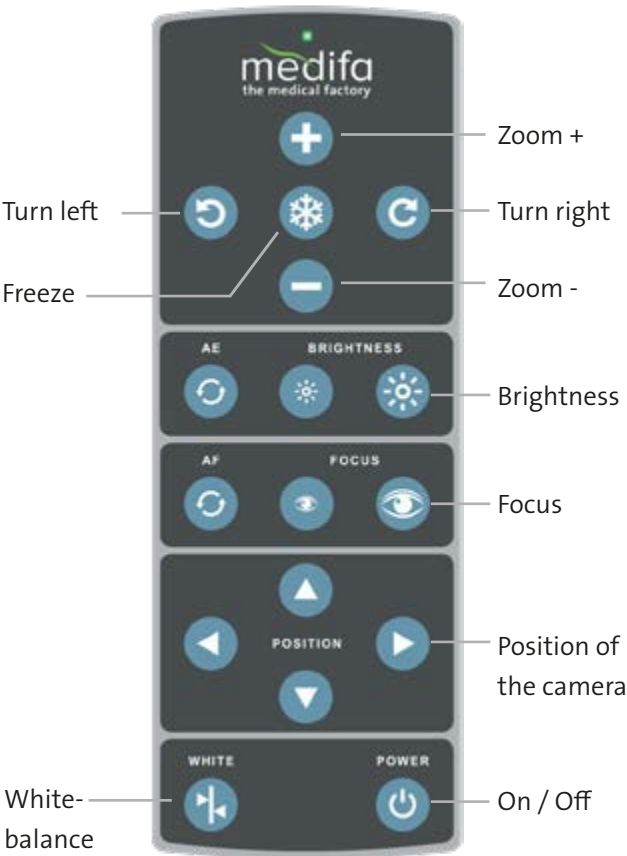

H!Light for your operating theatre



Everything from one source for your OR:

1. Operating lights
2. Modular room system RooSy
3. Mobile operating tables and accessories
4. Medical furniture and customized solutions

Technical data operating light	
Light intensity in Lux	40,000 lx - 160,000 lx
Adjustability of the light intensity	in 6 steps
Colour temperature (Kelvin)	3,000 - 5,500
Colour rendering index R _a	96
Colour rendering index R _g (blood)	95
Colour rendering index R ₁₃ (skin)	98
Focusable light field size (Diameter d ₁₀ in cm)	18 - 32 cm
Adjustability of the light field size	in 6 steps
d ₅₀ /d ₁₀ ratio	0.67
Working distance (in cm)	70 - 150 cm
Diameter of light head (in cm)	60 cm
Temperature increase in the head area	<0.5°
Number of LEDs	72
Life-span of the LEDs	> 60,000h
Depth of light field	65
L1+L2 at 60% Ec(cm)	
Depth of light field	120
L1+L2 at 20% Ec(cm)	
Total power consumption	160 VA

All values are conform to IEC 60601-2-41 and are subject to a general industrial tolerance of ± 10%.
*In accordance with the IEC 60601-2-41 standard, the maximum illuminance must not exceed 160,000 lx. The industrial tolerance is therefore only -10%.

Technical data wireless-camera	wireless HD	wireless 4K
	(Art.-no. 103200)	(Art.-no. 103200)
Resolution pixel	ca. 2,140,000 Pixel	8,510,000 Pixel
Resolution	1,920 x 1,080	3,840 x 2,160
Zoom	30x optical, 12x digital	20x optical, 12x digital
Image rotation	limitless, motorised	limitless, motorised
Video outputs	HDMI 1.4, DVI	HDMI 1.4, DVI
Aspect ratio	16:9	16:9
Dimensions (Ø, length)	68 mm x 179 mm	68 mm x 179 mm
Weight camera	1 kg	1 kg
Security of data transmission	exclusively between coupled transmitter + receiver AES 128 bit encryption	exclusively between coupled transmitter + receiver AES 128 bit encryption




Use our wireless camera systems in the operating theatre for detailed, precise representation and highest image quality

Our cameras offer you finely structured images, optimum detail and brilliant reproduction with clear, sharp contours and perfect colour reproduction from every position.

A wired solution is also available on request.

Configuration examples


Triple combination with monitor holder

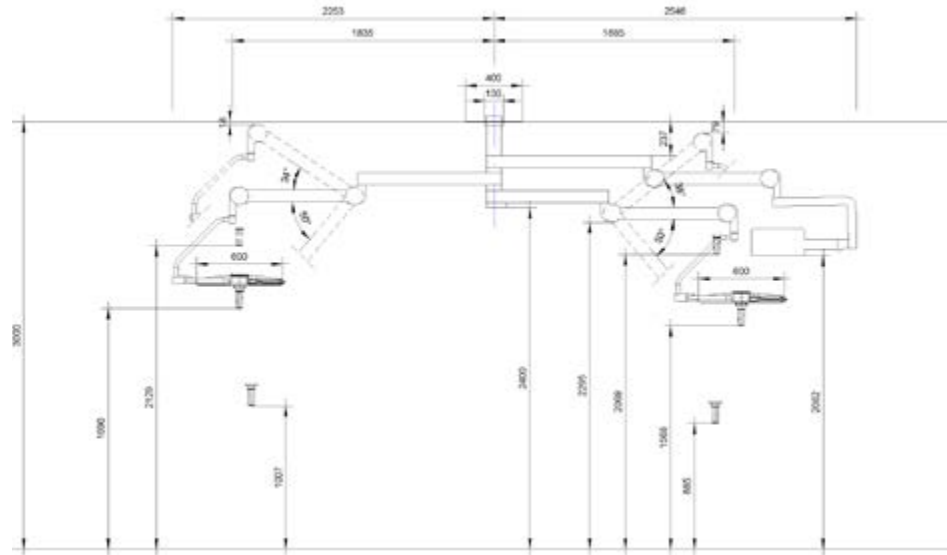

Combination of

- **Art.-no. 105075** Flat cover 400x400
- **Art.-no. 105060** Ceiling tube 105 mm
- **Art.-no. 105300** Central axis trio
- **Art.-no. 100200** medifa H!Light
- **Art.-no. 100400** medifa H!Light plus
- **Art.-no. 105700** 2x medifa handle adapter standard or
Art.-no. 105750 medifa handle adapter sterile control
- **Art.-no. 102400-1** medifa H!Light monitor 4K up to 32''
- **Art.-no. 102200-1** medifa H!Light monitor stop

Double combination

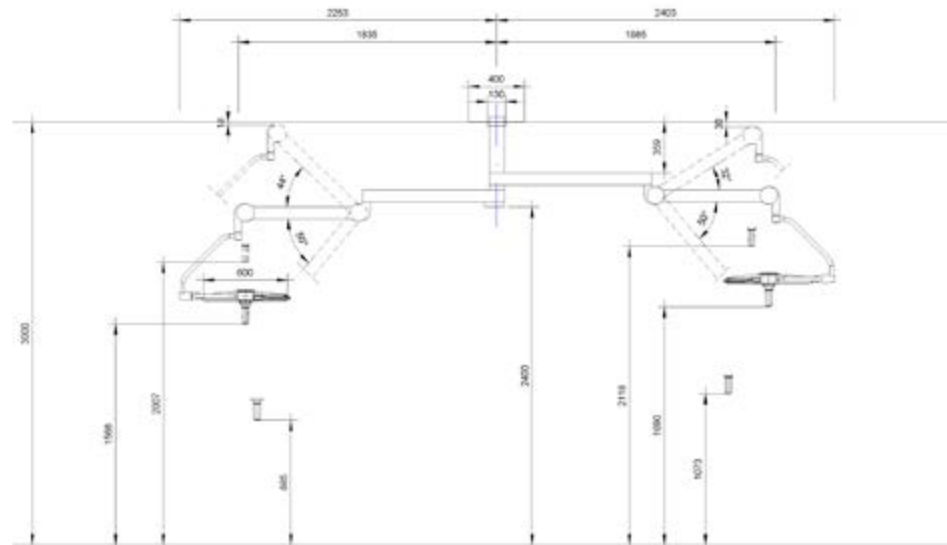

Combination of

- **Art.-no. 105075** Flat cover 400x400
- **Art.-no. 105060** Ceiling tube 105 mm
- **Art.-no. 105200** Central axis duo
- **Art.-no. 100200** medifa H!Light
- **Art.-no. 100400** medifa H!Light plus
- **Art.-no. 105700** 2x medifa handle adapter standard or
Art.-no. 105750 medifa handle adapter sterile control

Accessories



Ceiling fixation set Light

Art.-no. 105050


- For intermediate ceiling
- Max. distance intermediate/bare ceiling 1,200 mm for installation with flat cover
- Consisting of: 1x Ceiling anchor plate M16, 6x Distance unit M16



Ceiling tube 105 mm

Art.-no. 105060


- Ceiling tube diameter: 105 mm
- Standard length: 300 mm or 500 mm
- Other length on request
- Fixation hole pattern 6x M16 270 mm

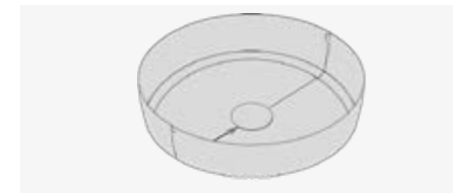

Cover 590

Art.-no. 105095

- Solo cover round
- D = 590 mm H = 125 mm
- Two part with 6 pcs. Mounting screws
- one part holding ring

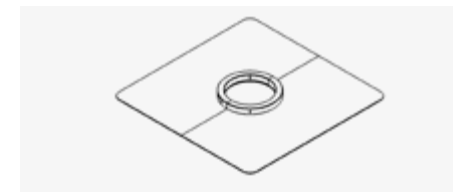

Flat cover 400x400

Art.-no. 105075

- Flat cover Solo
- 400 x 400 mm H = 2 mm
- Two part holding ring




Installationkit Laminar Air Flow

Art.-no. 105080

- Two units, two part holding ring


Accessories

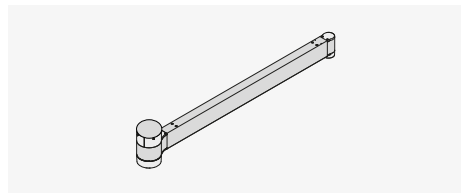

Central axis single

Art.-no. 105100


- Central axis and extension arm
- Extension arm length standard: 1,075 mm

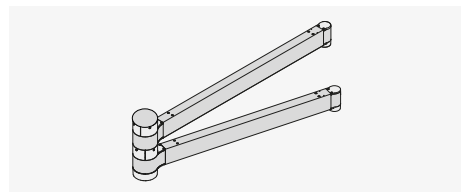

Central axis duo

Art.-no. 105200


- Central axis and extension arm
- Extension arm length standard: 925 mm / 1,075 mm

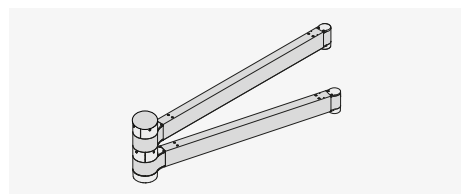

Central axis duo hybrid

Art.-no. 105210


- Central axis and extension arm
- Extension arm length standard: 1,375mm / 1,525 mm

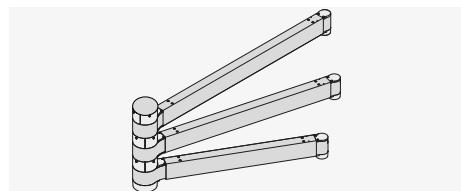

Central axis trio

Art.-no. 105300

- Central axis and extension arm
- Extension arm length standard: 775 mm / 925 mm / 1,075 mm

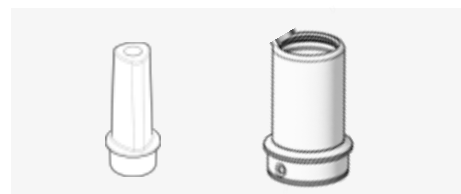

medifa handle standard (set of 5)

Art.-no. 104100

- Sterilisable handle

medifa handle camera (set of 5)

Art.-no. 104200

- Sterilisable handle camera

medifa

The medifa healthcare group unites competence concerning operating rooms under one roof.

medifa is an owner-run and internationally active complete solution provider in the medical technology.

The product and service portfolio covers all operating room requirements from modular wall, door, ceiling and cabinet systems (RooSy® Room Systems) for the turnkey operating room, as well as mobile operating tables, operating table accessories and operating lights, right up to examination chairs and medical

furniture. In addition, one of the largest German OEM suppliers in the medical technology sector also belongs to medifa.

medifa GmbH & Co. KG (Finnentrop), medifa hygienic rooms GmbH (Ötigheim) and medifa metall und medizintechnik GmbH in Rastatt are all part of this group of companies. The medifa healthcare group employs approximately 350 employees at three locations within Germany. The planning, engineering and production


take place exclusively in Germany, and with over 90% added value in the group of companies, the group guarantees fast reaction times.

The slogan “we care” means for us that we consistently orientate ourselves according to the requirements of our customers. Customers all over the world can rely on top-class consultation, processing and quality from a single source, made in Germany according to certified national and international quality standards.

Products

- Modular operating room systems: walls, doors, ceilings and cabinets (RooSy)
- Mobile operating tables and accessories
- Operating lights
- Examination and treatment chairs
- Medical furniture and trolleys
- Customized medical solutions

medifa worldwide

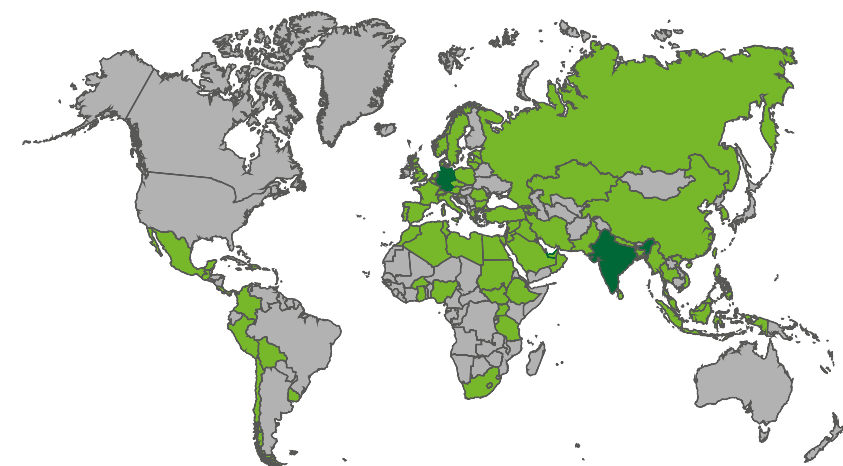

we care.



Version 05-20-EN

Subject to modifications due to technical development.
This document is intended to provide information to an
international audience outside of the US.

www.medifa.com

we care. Code



medifa GmbH & Co. KG

Industriestraße 5
57413 Finnentrop
Germany
fon +49 2721 71 77 -0
fax +49 2721 71 77 -255
info@medifa.com

we care.