

## Quick. Clever. Clip-on. Easy-to-attach infusion ports from A.M.I.

A.M.I.'s revolutionary Clip-a-Port System offers a variety of ports to meet every need, all featuring an integrated clip on the base of the port for fast, easy attachment to tissue.

- Two different catheter sizes
- Standard or slimline design for higher patient comfort

# A.M.I.® Clip-a-Port

The A.M.I. Clip-a-Port System is based on a unique design of clip-on infusion port for the venous supply of pharmaceuticals, fluid and nutrient solutions.

## Clip-on port


All models have a specially-designed base plate with an integrated clip, which can be slid under the fascia to fasten the port - no sutures and no time-wasting! A.M.I. Clip-a-Port Systems are also suitable for power injection at a maximum recommended flow rate of 5 ml/s and pressure not exceeding 325 psi (2240 kPa).

- Large silicone septum for easier needle insertion
- Dual safety edges prevent needle slippage
- Radiopaque catheter with length markings
- Very soft catheter tip
- Bio-compatible materials (titanium, silicone and poly-urethane)
- MR-compatible up to 3 Tesla; latex-free

## Port for venous cutdown

- Port with catheter protection
- Single-lumen catheter of polyurethane
- Vein hook
- Tunneller
- Rinsing needle
- A.M.I. Port Needle
- 5 cc syringe

7-piece set, delivered sterile

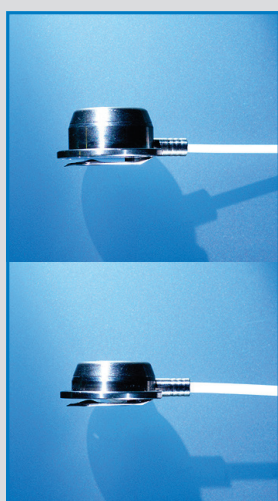

### Technical details

Diameter of silicone membrane:	12 mm
Length of catheter:	500 mm

### Standard design

Height of port (without clip):	13 mm
Catheter size F6:	CPS6751
Catheter size F8:	CPS6791

### Slimline design

Height of port (without clip):	9.8 mm
Catheter size F6:	CPS6771
Catheter size F8:	CPS6711





# A.M.I.® Clip-a-Port



## Port for venous cutdown or Seldinger technique

- Port with catheter protection
- Single-lumen catheter of polyurethane
- Vein hook
- Tunneller
- Rinsing needle
- A.M.I. Port Needle
- 5 cc syringe
- Seldinger puncture needle 18 G
- Dilator / peel-away sheath for F6 or F8 catheters
- Guidewire

10-piece set, delivered sterile

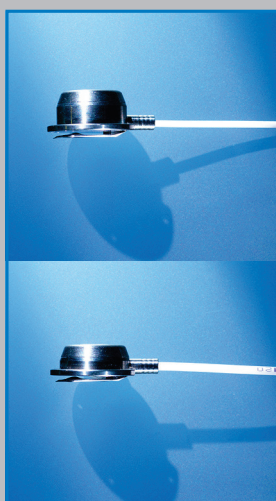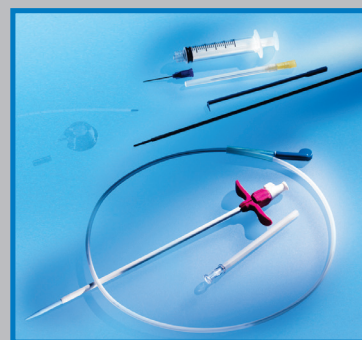

### Technical details

Diameter of silicone membrane: 12 mm  
Length of catheter: 500 mm

### Standard design

Height of port (without clip): 13 mm  
Catheter size F6: CPS6741  
Catheter size F8: CPS6781

### Slimline Design

Height of port (without clip): 9.8 mm  
Catheter size F6: CPS6761  
Catheter size F8: CPS6701

# A.M.I.® Clip-a-Port

Order Code	Product	Technical Details
CPS6751	<b>Clip-a-Port F6</b> - Standard-height infusion port of titanium with integrated clip and F6 catheter	1 set of 7 pieces (Cutdown) Delivered sterile
CPS6791	<b>Clip-a-Port F8</b> - Standard-height infusion port of titanium with integrated clip and F8 catheter	1 set of 7 pieces (Cutdown) Delivered sterile
CPS6771	<b>Clip-a-Port slimline F6</b> - Low-profile infusion port of titanium with integrated clip and F6 catheter	1 set of 7 pieces (Cutdown) Delivered sterile
CPS6711	<b>Clip-a-Port slimline F8</b> - Low-profile infusion port of titanium with integrated clip and F8 catheter	1 set of 7 pieces (Cutdown) Delivered sterile
CPS6741	<b>Clip-a-Port F6 kit</b> - Standard-height infusion port of titanium with integrated clip and F6 catheter - Catheter F6 Introducer Set	1 set of 10 pieces (Cutdown or Seldinger) Delivered sterile
CPS6781	<b>Clip-a-Port F8 kit</b> - Standard-height infusion port of titanium with integrated clip and F8 catheter - Catheter F8 Introducer Set	1 set of 10 pieces (Cutdown or Seldinger) Delivered sterile
CPS6761	<b>Clip-a-Port slimline F6 kit</b> - Low-profile infusion port of titanium with integrated clip and F6 catheter - Catheter F6 Introducer Set	1 set of 10 pieces (Cutdown or Seldinger) Delivered sterile
CPS6701	<b>Clip-a-Port slimline F8 kit</b> - Low-profile infusion port of titanium with integrated clip and F8 catheter - Catheter F8 Introducer Set	1 set of 10 pieces (Cutdown or Seldinger) Delivered sterile
CPS6941	<b>Catheter F6 Introducer Set</b> Consists of: - Seldinger puncture needle 18 G - Dilator / peel-away sheath for F6 - Guidewire	1 set Delivered sterile
CPS6951	<b>Catheter F8 Introducer Set</b> Consists of: - Seldinger puncture needle 18 G - Dilator / peel-away sheath for F8 - Guidewire	1 set Delivered sterile

## A.M.I. Headquarters:

A.M.I. GmbH  
Im Letten 1  
6800 Feldkirch . Austria  
t +43 5522 90505-0  
f +43 5522 90505-4006  
e info@ami.at  
[www.ami.at](http://www.ami.at)

