




IRRIGATION SET
ASPIRATION SET



VECTEC®
SINGLE USE
LAPAROSCOPIC
DEVICES



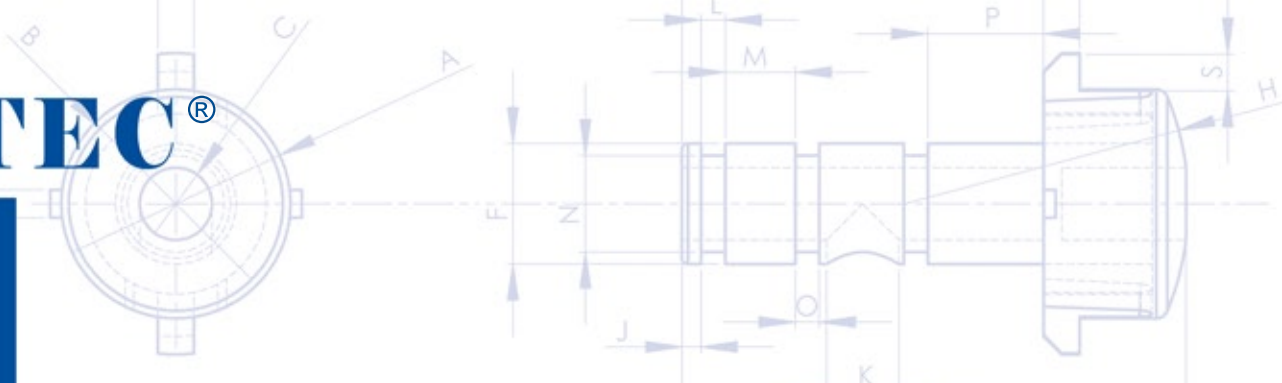

GRAVITY BASED IRRIGATION/ASPIRATION SET to fit general and laparoscopic surgery

PRODUCTS

Reference	Designation	Diameter	Length	Packaging
CAN250 5	irrigation/aspiration set	5 mm	330 mm	10 units per box
CAN250 5P	irrigation/aspiration set with anti-drip tip	5 mm	330 mm	10 units per box
CAN250 5W	irrigation/aspiration set with adaptor	5 mm	330 mm	10 units per box
CAN250 5GL	irrigation/aspiration set	5 mm	450 mm	10 units per box
CAN250 5DP	irrigation/aspiration set with double spike	5 mm	330 mm	10 units per box
CAN250 10	irrigation/aspiration set	10 mm	330 mm	10 units per box



Spike connector



Adaptor



Ø5 mm/10 mm cannula

ACCESSORIES

Reference	Designation	Diameter	Length	Packaging
CAN 10	Cannula only	10 mm	330 mm	10 units per box

The single use suction irrigation set consists in :

- An atraumatic interchangeable 5 mm or 10 mm sheath with distal lateral holes.
- An ergonomic handle fitted with double trumpet-valves with colour coded plunger, also suitable for biopsy
- A transparent and flexible 2,5 meter long double tubing set
- A single spike connection for irrigation bags and an on-and-off clamp
- A soft universal connector fitting to any suction canister

MD Class IIa
Sterilized by ETO

See instructions for use
Manufactured by Promepla

CE
0459

G
MED
NF EN ISO 13485



Distributed by:
Péters Surgical
42, rue Benoit Frachon
93013 Bobigny Cedex - FRANCE
Tel: +33 (0)1 48 10 62 62
Fax: +33 (0)1 48 10 22 99
commercial@peters-surgical.com
www.peters-surgical.com