



# Edwards Fogarty<sup>®</sup> Catheters


The only name in clot  
management backed  
by over 50 years of  
worldwide experience.



Edwards

# For over 50 years, clinicians have come to address the real problem of clot

The reason is simple. Through collaboration and clinical learning, Edwards Lifesciences has continued to evolve and advance its Fogarty® catheters to meet the exacting needs of specialists for vascular procedures. Their proven performance and manufacturing quality has led surgeons to specify genuine Fogarty® catheters – available exclusively from Edwards – for millions of procedures worldwide.

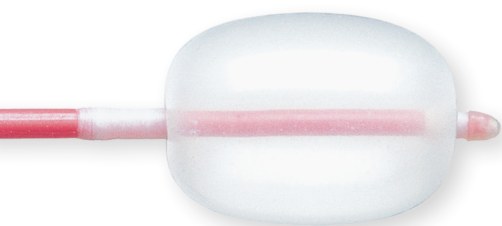

## Fresh, soft emboli/thrombi

### **Fogarty® arterial embolectomy catheter**

Removal of fresh, soft emboli and thrombi from vessels in the arteries of the non-central circulatory system


- Hand-tied using a recessed winding technique to secure balloon to catheter
- Balloon exhibits a symmetry that exerts uniform contact with vessel walls for even pressure and precise traction
- Rounded tip promotes easy insertion for reduced trauma



### **Fogarty® Fortis arterial embolectomy catheter**

Removal of fresh, soft emboli and thrombi from vessels in the arteries of the non-central circulatory system

- Can withstand greater pull-force than the Fogarty® arterial embolectomy catheter
- The addition of barium to the catheter tubing provides radiopacity

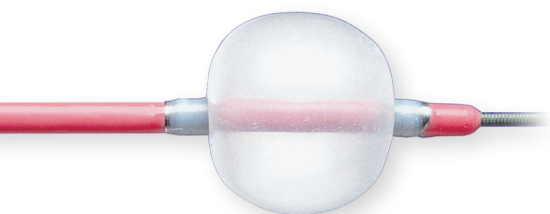

## Fresh, soft emboli/thrombi, temporary occlusion, infusion

### **Over-the-wire Fogarty® thru-lumen embolectomy catheter**

Multi-use guidewire-compatible design enables removal of fresh, soft emboli and thrombi from vessels in the arterial system as well as temporary occlusion of blood vessels, infusion of fluids, and blood sampling

- Fogarty® balloon conforms to the vessel to maximize clot removal while minimizing trauma
- Provides the solution for distal thrombus removal in the arterial system
- The addition of barium to the catheter shaft material and stainless steel markers proximal and distal to the balloon provide radiopacity\*

\*3F Thru-Lumen embolectomy catheter body does not include barium

# nted on real Fogarty® catheters management.

## Mature, adherent thrombi

### Fogarty® Corkscrew catheter

Indicated for removal of emboli and thrombi from either native vessels or synthetic grafts in the arterial system

- For extracting thrombus that is too resistant for an elastomeric balloon
- Corkscrew design increases the surface area for entrapping fibrous material
- Radiopaque and flexible

### Fogarty® graft thrombectomy catheter

Designed for removing organized thrombus and adherent thrombotic material from synthetic bypass grafts, specifically, AV fistulas for dialysis access and grafts for iliofemoral bypass

- For thrombus that is too resistant for an elastomeric balloon
- Provides a larger diameter range than the Fogarty® Corkscrew catheter
- The pulling force can be varied by adjusting the diameter

## Venous thrombi

### Fogarty® venous thrombectomy catheter

Designed for removing thrombi in the venous system

- The long, soft tip permits advancement past the venous valves without undue trauma



## You'll want a real

A close-up of the tip of a Fogarty catheter, showing the brand name "Fogarty®" printed on the white, conical tip.

**Fogarty®**

## Fogarty® clot management catheters

Model	Length (cm)	Catheter size (F)	Membrane/spiral diameter max. extended (mm)	Membrane/spiral diameter max. retracted (mm)
<b>Fogarty® Corkscrew catheters</b>				
140806	80	4	2.7	6
140808	80	5	3.1	8
1408010	80	6	3.3	10
<b>Fogarty® graft thrombectomy catheters</b>				
160245F	50	5	5.0	16
160246F	50	6	6.0	18

Model	Length (cm)	Catheter size (F)	Inflated balloon diameter (mm)	Max. size of deflated balloon (F)	Max. liquid capacity (ml)	Max. recommended pull force (lbs.)	Compatible introducer (F)	Compatible guidewire (in.)
<b>Fogarty® arterial embolectomy catheters (pouches, 5 pack)</b>								
120602F	60	2	4	3.4	–	0.5	3.5	–
120403F	40	3	5	4.3	0.20	0.7	4	–
120803F	80	3	5	4.3	0.20	0.7	4	–
120404F	40	4	9	5.0	0.75	1.5	5	–
120804F	80	4	9	5.0	0.75	1.5	5	–
120805F	80	5	11	6.0	1.50	2.0	6	–
120806F	80	6	13	7.0	2.00	2.5	7	–
120807F	80	7	14	8.0	2.50	3.5	8	–
<b>Fogarty® Fortis arterial embolectomy catheters (pouches, 5 pack)</b>								
120804FF	80	4	9	5.0	0.75	–	5	–
<b>Over-the-wire Fogarty® thru-lumen embolectomy catheters</b>								
12TLW803F	80	3	5	4.5	0.15	0.7	4.5	0.018
12TLW404F	40	4	9	5.0	0.50	1.5	5	0.025
12TLW804F	80	4	9	5.0	0.50	1.5	5	0.025
12TLW805F35	80	5.5	11	6.0	0.90	2.0	6	0.035
12TLW406F	40	6	13	7.0	1.40	2.5	7	0.035
12TLW806F	80	6	13	7.0	1.40	2.5	7	0.035
12TLW807F	80	7	14	8.0	1.60	3.5	8	0.038
<b>Fogarty® venous thrombectomy catheters</b>								
320806F	80	6	12	7.7	1.50	–	–	–
320808F	80	8	13	8.8	2.25	–	–	–
32080810F	80	8	19	10.0	4.00	–	–	–

Genuine Fogarty® catheters are available exclusively from Edwards Lifesciences. To learn more, visit [edwards.com/fogarty](https://www.edwards.com/fogarty).

For professional use. For a listing of indications, contraindications, precautions, warnings, and potential adverse events, please refer to the Instructions for Use (consult [eifu.edwards.com](https://www.edwards.com) where applicable).

Edwards devices placed on the European market meeting the essential requirements referred to in Article 3 of the Medical Device Directive 93/42/EEC bear the CE marking of conformity.

Edwards, Edwards Lifesciences, the stylized E logo, Fogarty, Fogarty Corkscrew, and Fogarty Fortis are trademarks of Edwards Lifesciences Corporation or its affiliates. All other trademarks are the property of their respective owners.

© 2020 Edwards Lifesciences Corporation. All rights reserved. PP--EU-0084 v1.0

Edwards Lifesciences • Route de l'Etraz 70, 1260 Nyon, Switzerland • [edwards.com](https://www.edwards.com)



Edwards